

Validation of the Pulmonary Function System for Use on the International Space Station

Frank A. McCleary, M.S
Wyle Life Sciences, Houston, TX, USA

Alan D. Moore, Jr., Ph.D.
Wyle Life Sciences, Houston, TX, USA

R. Donald Hagan, Ph.D.
NASA Johnson Space Center, Houston, TX, USA

THE NASA STI PROGRAM OFFICE . . . IN PROFILE

Since its founding, NASA has been dedicated to the advancement of aeronautics and space science. The NASA Scientific and Technical Information (STI) Program Office plays a key part in helping NASA maintain this important role.

The NASA STI Program Office is operated by Langley Research Center, the lead center for NASA's scientific and technical information. The NASA STI Program Office provides access to the NASA STI Database, the largest collection of aeronautical and space science STI in the world. The Program Office is also NASA's institutional mechanism for disseminating the results of its research and development activities. These results are published by NASA in the NASA STI Report Series, which includes the following report types:


- **TECHNICAL PUBLICATION.** Reports of completed research or a major significant phase of research that present the results of NASA programs and include extensive data or theoretical analysis. Includes compilations of significant scientific and technical data and information deemed to be of continuing reference value. NASA's counterpart of peer-reviewed formal professional papers but has less stringent limitations on manuscript length and extent of graphic presentations.
- **TECHNICAL MEMORANDUM.** Scientific and technical findings that are preliminary or of specialized interest, e.g., quick release reports, working papers, and bibliographies that contain minimal annotation. Does not contain extensive analysis.
- **CONTRACTOR REPORT.** Scientific and technical findings by NASA-sponsored contractors and grantees.

- **CONFERENCE PUBLICATION.** Collected papers from scientific and technical conferences, symposia, seminars, or other meetings sponsored or cosponsored by NASA.
- **SPECIAL PUBLICATION.** Scientific, technical, or historical information from NASA programs, projects, and mission, often concerned with subjects having substantial public interest.
- **TECHNICAL TRANSLATION.** English-language translations of foreign scientific and technical material pertinent to NASA's mission.

Specialized services that complement the STI Program Office's diverse offerings include creating custom thesauri, building customized databases, organizing and publishing research results . . . even providing videos.

For more information about the NASA STI Program Office, see the following:

- Access the NASA STI Program Home Page at <http://www.sti.nasa.gov>
- E-mail your question via the internet to help@sti.nasa.gov
- Fax your question to the NASA Access Help Desk at (301) 621-0134
- Telephone the NASA Access Help Desk at (301) 621-0390
- Write to:
NASA Access Help Desk
NASA Center for AeroSpace Information
7115 Standard
Hanover, MD 21076-1320

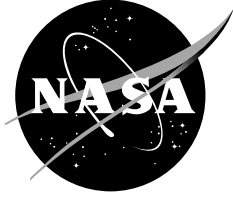

Validation of the Pulmonary Function System for Use on the International Space Station

*Frank A. McCleary, M.S.
Wyle Life Sciences, Houston, TX, USA*

*Alan D. Moore, Jr., Ph.D.
Wyle Life Sciences, Houston, TX, USA*

*R. Donald Hagan, Ph.D.
NASA Johnson Space Center, Houston, TX, USA*

Available from:

NASA Center for Aerospace Information
7121 Standard Drive
Hanover, MD 21076-1320
301-621-0390

National Technical Information Service
5285 Port Royal Road
Springfield, VA 22161
703-605-6000

This report is also available in electronic form at <http://ston.jsc.nasa.gov/collections/TRS/>

TABLE OF CONTENTS

1.0	Introduction.....	1
1.1	Purpose	2
2.0	Methods.....	2
2.1	Subjects.....	2
2.2	Protocol.....	3
2.3	Metabolic Gas Analysis Systems Description	4
2.4	Metabolic Gas Analysis Data Reduction.....	4
2.5	Data Analysis.....	5
3.0	Results	5
3.1	Oxygen Consumption (VO_2)	5
3.2	Remaining Metabolic Gas Analysis Measures	7
4.0	Discussion.....	8
5.0	References	10

LIST OF TABLES

Table 1. Subject Characteristics.....	2
---------------------------------------	---

LIST OF FIGURES

Figure 1. Peak Cycle Protocol.....	3
Figure 2. Comparison of ParvoMedics reference system and PFS measured VO_2 at the work rates that all subjects could complete.	6
Figure 3. Relationship between all exercise VO_2 measurements collected by the ParvoMedics reference system and those made by the PFS.	6
Figure 4. Comparison of ParvoMedics and PFS measured VCO_2 values.....	7
Figure 5. Comparison of ParvoMedics and PFS measured expired pulmonary ventilation values.	7
Figure 6. Comparison of ParvoMedics and PFS measured FEO_2 values.....	7
Figure 7. Comparison of ParvoMedics and PFS measured $FECO_2$ data.	7

1.0 INTRODUCTION

Aerobic deconditioning occurs during long duration spaceflight despite the use of exercise countermeasures (Convertino, 1996). As a part of International Space Station (ISS) medical operations, periodic tests designed to estimate aerobic capacity are performed prior to, during, and after missions of greater than 30 days in duration. These tests track changes in aerobic fitness and determine the effectiveness of exercise countermeasures.

Crew members selected for missions aboard the ISS perform a graded exercise test on a cycle ergometer approximately 270 days prior to their scheduled launch date in order to measure peak oxygen consumption (VO_{2PK}) and peak heart rate (HR_{PK}). Approximately 30 to 45 days prior to launch, crew members perform a submaximal cycle ergometer test at work rates set to elicit 25, 50 and 75% of their preflight VO_{2PK} . This test, known as the Periodic Fitness Evaluation (PFE), serves as a baseline measure to which subsequent in- and post-flight exercise tests are compared. While onboard the ISS, crew members normally perform the PFE beginning with flight day (FD) 14 and every 30 days thereafter. The PFE is also conducted 5 and 30 days following flight.

Using PFE data, aerobic fitness is estimated by quantifying the oxygen consumption (VO_2) vs. heart rate (HR) relationship using linear regression and calculating the VO_2 that would occur at the crew member's previously measured HR_{PK} . Currently, for data collected during flight, this technique assumes that the pre- vs. in-flight oxygen consumption per given cycle workload is similar. However, the validity of this assumption is based upon a sparse amount of data collected during the Skylab era (Michel, et al., 1977). The method of using HR and cycle ergometer work rates has been used to estimate aerobic fitness in normal gravity (Astrand and Ryhming, 1954; Lee, 1993). Due to spaceflight-induced physiological alterations, such as shifts in extracellular fluid (e.g., plasma) volume, this method may not be valid during spaceflight. In addition, the ergometer onboard the ISS is vibration-isolated and moves with the astronaut's application of force into the pedals. The effect of this movement on the VO_2 of cycle exercise on ISS has not been quantified.

Though the measurement of VO_2 during ISS flight has not been conducted to date, it is a long-waived requirement found within the International Space Station Medical Operations Requirement Document (ISS MORD), Rev C. An attempt to meet this requirement was made using accessories for the Gas Analyzer System for Metabolic Analysis Physiology (GASMAP), a mass spectrometer developed by the NASA's ISS Human Research Facility (HRF) to support multiple studies onboard the ISS, including measurements of VO_2 during exercise. However, ground-based laboratory testing showed that, though the GASMAP performed well in the analysis of static gas samples, the software used to measure VO_2 was not accurate (A. Moore, personal communication, February 2000). Subroutines within the software were unable to calibrate the turbine flowmeter using a standard three-liter calibration syringe. Furthermore, the delay calculation required to measure breath-by-breath VO_2 was continuously out of phase and unable to align the expired gas fractions with the corresponding expired ventilation measurement. When compared to the Johnson Space Center Exercise Physiology Lab's (EXL) reference metabolic gas analysis system, the GASMAP was 24 and 11% higher for submaximal workloads of 25 and 50% VO_2 peak, respectively. It was concluded that without further design and software modifications, the GASMAP was not an adequate device for measuring VO_2 onboard the ISS. The expense of engineering support and hardware modifications required to configure the GASMAP to perform accurate measures of VO_2 was considered prohibitive at the time and were not implemented.

A contractor to the European Space Agency (ESA) recently developed a device capable of measuring exercise VO_2 onboard the ISS. The device, known as the Pulmonary Function Module/Photoacoustic

Analyzer Module (PFM/PAM), is part of the ISS Pulmonary Function System (PFS). For the purposes of this report, the PFM/PAM and PFS are referred to simply as the PFS. The PFS is located in the Human Research Facility (HRF) Rack 2 of the Destiny Laboratory onboard the ISS. The manufacturer of the PFS, Damec (Damec Research ApS, Odense, Denmark), has a history of developing equipment designed to measure the effects of microgravity on the human respiratory system. Previous products by this company include the Respiratory Monitoring System (RMS) for the Anthrorack used during the Space Lab missions, the RMS-II used during the EuroMir 95 mission, and the Advanced Respiratory Monitoring System used aboard STS-107.

The PFS was initially launched to the ISS aboard STS-114 (LF 1) in July 2005, followed by a hardware upgrade that was launched on a Russian Progress supply vehicle in April 2006. A collaborative effort has been initiated between NASA and ESA Life Sciences personnel to integrate the use of the PFS into the PFE tests for the support of Medical Operations objectives. Jensen, et al. (2002) performed an in-depth comparison of the AMIS 2001, which is the commercial version of the PFS, versus the standard Douglas bag method and found that the system was reliable and accurate for measuring VO₂. However, prudence dictates that an evaluation of the flight-like PFS device be required to determine if it provides valid metabolic gas analysis prior to implementation onboard the ISS.

1.1 Purpose

The purpose of this investigation was to compare exercise metabolic gas analysis measurements (including VO₂) obtained by the PFS to those collected using a reference metabolic gas analysis system: the ParvoMedics TrueOne[®] 2400 system (ParvoMedics, Salt Lake City, UT). This system has been extensively validated (Basset, et al., 2001; Crouter, et al., 2006) and is currently utilized by the NASA's Exercise Physiology Laboratory for pre- and post-flight testing astronauts assigned to ISS flights.

2.0 METHODS

2.1 Subjects

Eight healthy male subjects volunteered to perform two peak cycle tests over a 12-day period (Table 1). All subjects completed a modified Air Force Class III physical exam prior to participation and received written and verbal explanations of test protocols before providing written informed consent. The NASA-Johnson Space Center Committee for the Protection of Human Subjects reviewed and approved the test protocols and procedures.

Table 1: Subject Characteristics (mean ± SD)

Characteristic	n = 8
Age (yrs)	33 ± 6
Height (cm)	181.9 ± 5.0
Weight (kg)	82.2 ± 9.1
VO _{2PK} (l·min ⁻¹)	4.55 ± 0.82
VO _{2PK} (ml·kg ⁻¹ ·min ⁻¹)	55.9 ± 14.9

2.2 Protocol

Each subject performed two peak cycle ergometer tests using a LODE Excalibur Sport (LODE, AN Groningen, The Netherlands). Metabolic gas analysis was accomplished using one of two systems in random order: the ParvoMedics TrueOne[®] 2400 system or the PFS. Subjects completed the two tests within a 12-day period. A minimum of 48 hours was allowed between each test to prevent residual soreness and fatigue. Sessions were conducted at approximately the same time of day, limiting variance to within ± 2 hours of the initial testing session to limit any circadian effect on physiological variance (Carter, 2002). Subjects were also requested to maintain the same dietary and sleeping habits throughout the testing period. Compliance was monitored through a subject screening form completed by the subject before each test.

Metabolic gas analysis was conducted continuously throughout the cycle test protocol (Figure 1). The cycle work rate was increased according to the protocol until the subject indicated that he could no longer continue. Each subject was allowed to choose his desired pedal cadence (revolutions per minute [RPM]) and was required to maintain that same cadence throughout the two tests. Lepers, et al. (2001) found no difference in VO_2 when subjects were allowed to choose their optimal cadence between 70 and 100 RPMs. All test subjects in the current evaluation maintained a pedal cadence within a range of 80 to 100 RPMs.

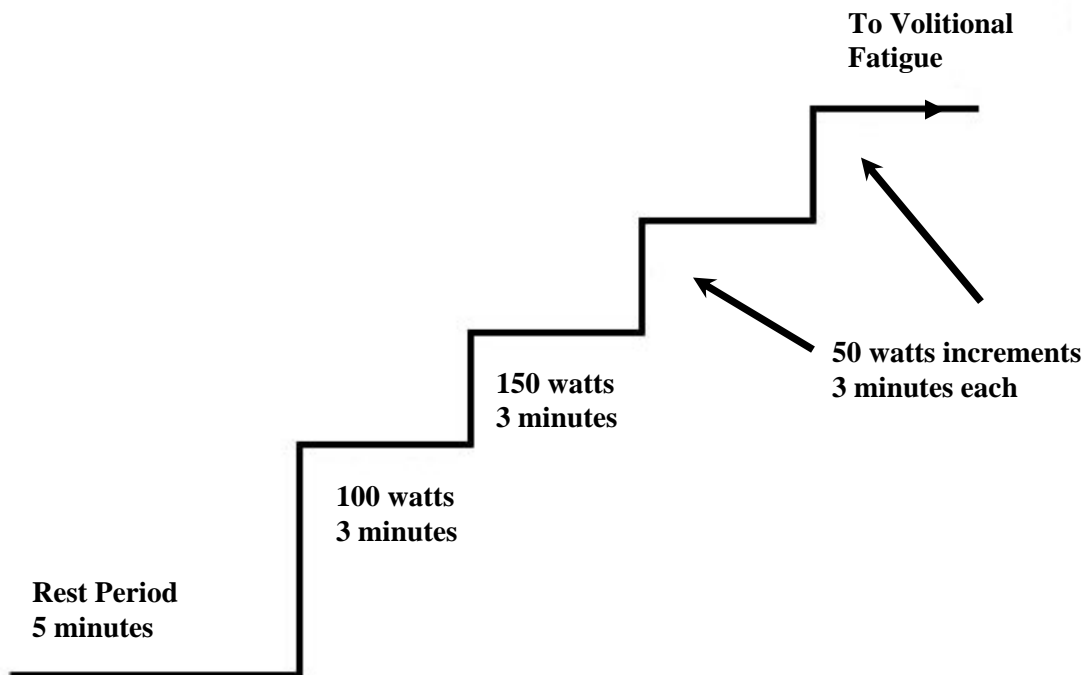


Figure 1: Peak Cycle Protocol

HR was measured electrocardiographically (Q-Stress, Quinton Instruments, Seattle, WA) and blood pressure was measured using a mercurial sphygmomanometer and stethoscope. Ratings of Perceived Exertion (RPE) were reported during the last 30 seconds of each stage. Blood pressure was not recorded above 250 watts (W) if the subject showed a normal blood pressure response to increases in workload for workloads below 250 W. Immediately post exercise, blood pressure was recorded and the subject was monitored for any adverse effects caused by maximal exertion.

2.3 Metabolic Gas Analysis Systems Description

The ParvoMedics TrueOne[®] 2400 system uses a paramagnetic oxygen analyzer (operating range 0-25% O₂) and an infrared single beam, single wavelength carbon dioxide analyzer (operating range 0-15% CO₂) to measure the composition of expired gasses. The subject inspires through a two-way non-rebreathing valve (Hans Rudolph Model 2700, Kansas City, MO) and expired air composition is analyzed in a 4-liter mixing chamber. The inspired gas composition is assumed to be standard atmospheric values (i.e., 20.93% O₂ and 0.03% CO₂). Expired ventilation is measured using a Hans Rudolph Model 3813 linear pneumotach (operating flow range 0-800 L/min). Computational software is provided with the system and runs with the Windows XP Professional operating system.

The Damec PFS uses two types of technology for gas analysis. For carbon dioxide, a photoacoustic method of gas analysis is utilized. In this technique, the gas sample is exposed to intermittent infrared light. The gas sample absorbs the light and the absorbed energy results in a rise in pressure by heating. The intermittent infrared light is divided into different pulsation frequencies and is filtered optically. Each optical filter allows only a specific wavelength of light to pass through and the wavelengths correspond to the infrared absorption spectra of the sample gasses. When the light source is removed, the gas cools down, resulting in a pressure fluctuation. By choosing a pulsation frequency in the audible range, the pressure fluctuation becomes an acoustic signal which is detected by a microphone. The audible sounds recorded by the microphone are analyzed and the amplitude of each signal is used to calculate the gas concentration. The operating range for carbon dioxide is from 0-12%. Gasses other than carbon dioxide can be detected utilizing the photoacoustic method, but an evaluation of these was beyond the scope of our investigation. For oxygen analysis, an Oxigraf[™] sensor is used in the PFS. The Oxigraf[™] technology is based on a laser diode absorption spectroscopy technique. The sample gas is exposed to a laser tuned to a wavelength of 760 nm (the peak of oxygen absorption). The laser signal is attenuated in proportion to the concentration of oxygen present in the sample. The operating range is from 0-100%.

When the PFS is used during exercise testing, the subject inspires through a Damec custom-designed two-way non-rebreathing valve and the expired gasses are sampled in a 15-liter anesthesia bag, which serves as a mixing reservoir. Inspired gas concentrations are measured on the inspired side of the respiratory circuit. This is done because inspired gas concentrations onboard the spacecraft typically deviate slightly from normal atmospheric values. Ventilation is measured on the inspired side of the non-rebreathing valve using a Damec custom-designed pneumotach (operating flow range 0-900 L/min).

2.4 Metabolic Gas Analysis Data Reduction

Data were collected continuously by the ParvoMedics system and the PFS. Data were averaged in 30-second intervals to the nearest whole breath. VO_{2PK} was accepted as the highest VO_2 attained for a single 30-second period. The dependent variables in this investigation were VO_2 , carbon dioxide production (VCO_2), pulmonary ventilation (V_E), and fractions of expired oxygen and carbon dioxide (F_{EO_2} and F_{ECO_2}). Custom data reduction templates were created using Microsoft Excel (Microsoft Corporation, Redmond, WA).

Data collected using the ParvoMedics were automatically reduced to 30-second intervals by the ParvoMedics software. The exported data were further reduced using a data reduction template and averaged over the last 30-seconds of each stage. The ParvoMedics sampled the QRS signal directly from the Q-Stress ECG system and calculated HR. These HR data were reduced in the same fashion as the metabolic gas analysis variables.

Damec provided a Ground Support Equipment (GSE) software package to calculate the dependent measures from the raw data signals measured by the PFS. Similar to the ParvoMedics data, all PFS data were reduced in 30-second intervals and exported to Excel. A data cable to collect HR data from the Q-Stress ECG system directly by the PFS was not available; therefore, HR data were entered into the PFS data reduction template from the Q-Stress ECG printout.

2.5 Data Analysis

Data analysis for each dependent variable was conducted as a series of preplanned comparisons (dependent T-tests with a Bonferroni adjustment made for multiple comparisons). This data reduction strategy was selected instead of the more traditional analysis of variance for two reasons: 1.) The dependent measurements (VO_2 , VCO_2 , and V_E) will increase with exercise intensity, with no need for a statistical test to corroborate a physiological certainty. 2.) The primary question of interest involves determining if the measured metabolic gas analysis variables differed between the devices at any stage of exercise. Testing of this question is much more straightforward using preplanned comparisons. A significant difference was deemed to occur if $p < 0.05$. A linear regression analysis comparing the VO_2 data collected from both devices across all exercise stages was also performed and confidence intervals developed for the slope and intercept of the relationship. Finally, a comparison of the $\text{VO}_{2\text{PK}}$ data attained from each of the devices was performed using a paired T-test.

3.0 RESULTS

The average work rate attained by the subjects was 319 ± 53 W. All subjects completed each exercise stage through 250 W. Therefore, the planned comparisons of the dependent variables were conducted at 100, 150, 200 and 250 W.

3.1 Oxygen Consumption (VO_2)

VO_2 did not differ between devices at 100, 200 and 250 W; however, a statistically significant difference between the devices was observed at the 150 W stage (Figure 2). At this stage, the percent difference in VO_2 between the two devices was 5.4%. When the VO_2 values measured on the ParvoMedics and the PFS were compared using linear regression, the Pearson product-moment coefficient of correlation was high ($r=0.987$), indicating that 97.5% of the variance observed in the measure was shared between the two devices (Figure 3). In addition, the intercept of the relationship was not significantly different than 0 and the slope of the relationship did not differ from 1.0. The peak oxygen consumption of the subjects did not differ between the two systems (ParvoMedics: 4.47 ± 0.73 L/min; PFS: 4.64 ± 0.77 L/min, a 3.8% difference).

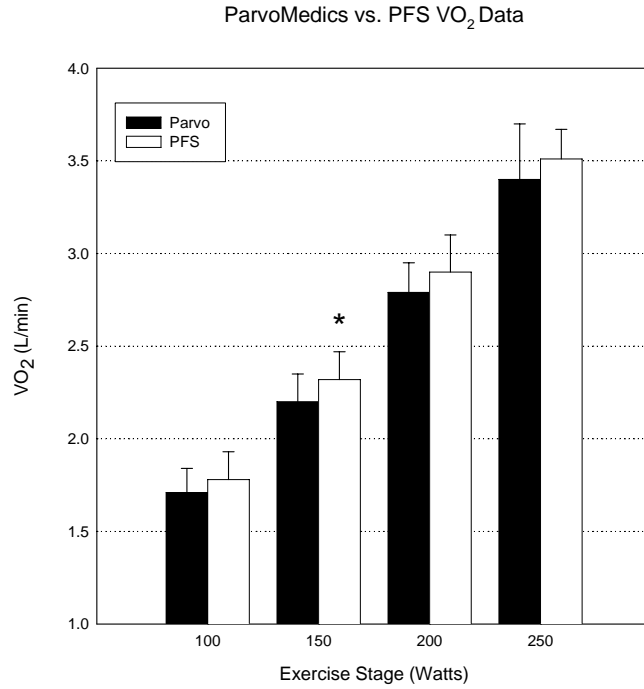


Figure 2: Comparison of ParvoMedics reference system and PFS measured VO₂ at the work rates that all subjects could complete. * (p<.05 – difference between devices at indicated stage)

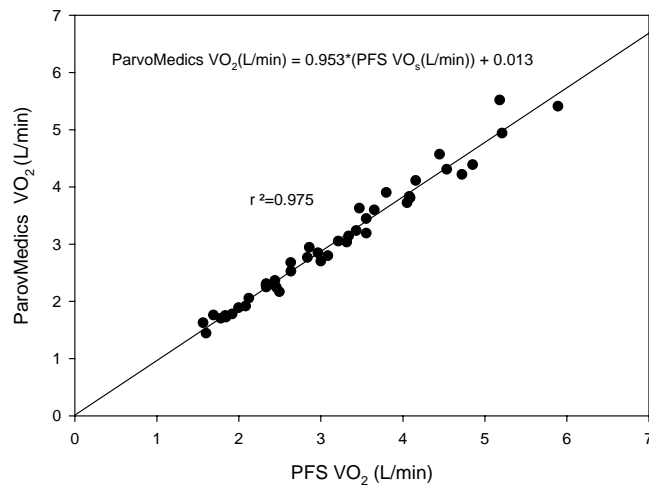


Figure 3: Relationship between all exercise VO₂ measurements collected by the ParvoMedics reference system and those made by the PFS.

3.2 Remaining Metabolic Gas Analysis Measures

None of the remaining metabolic gas analysis measures examined (i.e., VCO_2 , V_E , FEO_2 , $FECO_2$) demonstrated any statistically significant differences between the ParvoMedics reference system and the PFS within any exercise stage (Figures 4-7).

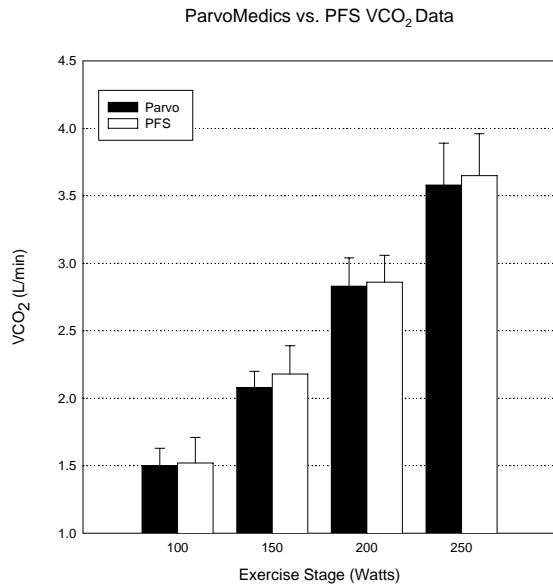


Figure 4: Comparison of ParvoMedics and PFS measured VCO_2 values.

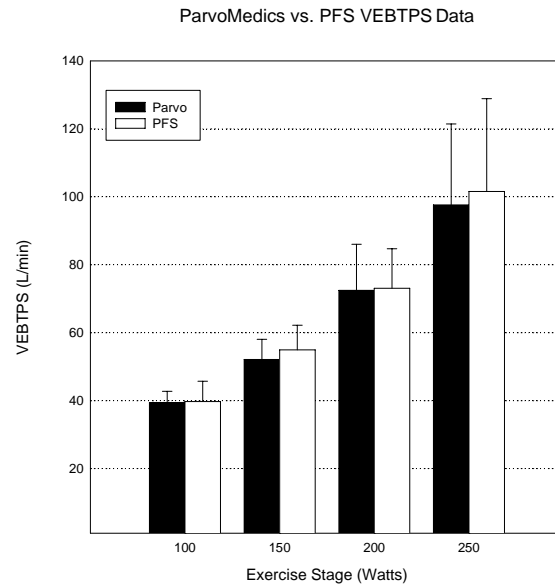


Figure 5: Comparison of ParvoMedics and PFS measured expired pulmonary ventilation values.

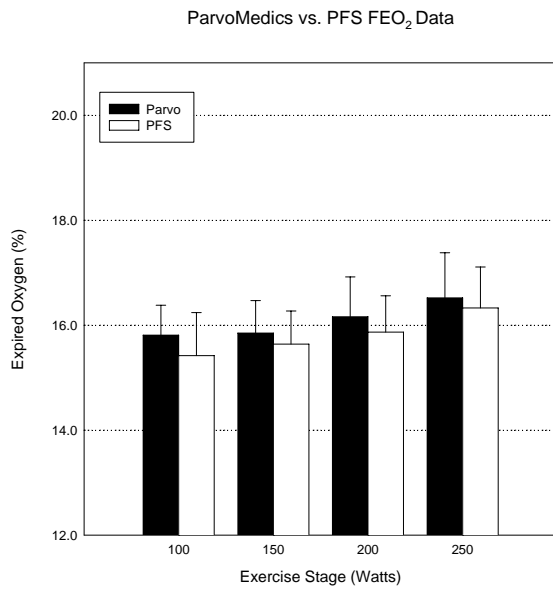


Figure 6: Comparison of ParvoMedics and PFS measured FEO_2 values.

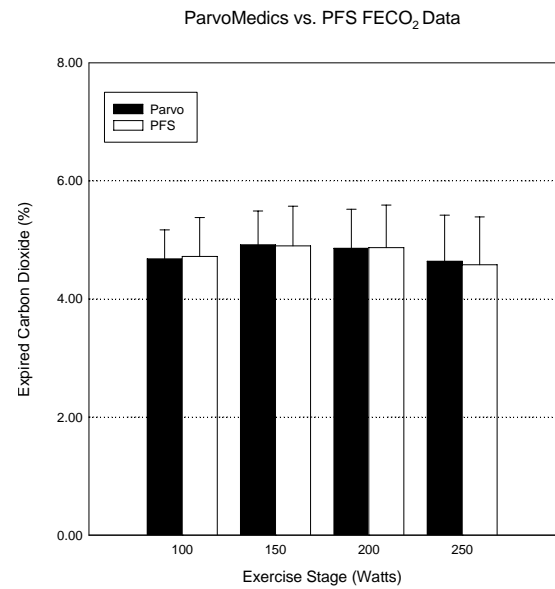


Figure 7: Comparison of ParvoMedics and PFS measured $FECO_2$ data.

4.0 DISCUSSION

Comparison of metabolic gas analysis systems is complicated by human variability in day-to-day exercise responses and performance. One way to control for variation in responses is to test both systems simultaneously, either in series or in parallel. While this is possible using some systems, the investigators of the current study and the designers of the hardware concluded that use of the ParvoMedics system and the PFS simultaneously may induce error related to a novel test setup. For example, if the ParvoMedics was sampling “in-line” past the PFS, the increased pressure of gas passing through the ParvoMedics may change the mixing bag characteristics of the PFS. Thus, two separate tests were performed on the subjects. While every attempt was made to limit variation, including limiting testing for subjects to the same time of day, it was reported that day-to-day variability in VO_2 ranges from 4 to 6% (Shephard, 1984), while others reported variability in $\text{VO}_{2\text{PK}}$ as high as 10 to 12% (Versteeg, et al., 1989).

The ParvoMedics reference system and the PFS measures of VO_2 did not statistically differ for the majority of exercise stages examined (Figure 2). For the 150 W stage, the systems statistically differed by 5.4%. This appears to be a random finding that has no apparent physiological or methodological explanation. One of the subjects did experience an unusually large difference between the two devices at that stage (2.49 L/min PFS vs. 2.17 L/min ParvoMedics, a 14.7% difference), and this undoubtedly contributed to the statistical difference. Consideration was given to removing the subject as an “outlier” from the data set, but because of the small sample size and the fact that all of the other stage data appeared normal for this subject, he was retained. From a clinical relevance standpoint, a 5.4% difference (i.e., the difference between the mean values of the two systems at 150 W), is within most of the day-to-day variations of VO_2 reported in the literature; therefore, based upon the stages of exercise examined, the two devices exhibited acceptable agreement. In addition, the peak VO_2 of the subjects demonstrated a non-significant difference of 3.8% between the two devices, also within the ranges reported in the literature of day-to-day variation in $\text{VO}_{2\text{PK}}$.

The VO_2 relationship between the two systems shows that they exhibit very close agreement. One would expect high agreement between two devices designed to measure the same physiological parameter, thus a high Pearson product-moment coefficient of correlation (r); however, a high r alone does not ensure that the results between the two devices are similar (Bland and Altman, 1986). If a high r is coupled with the slope of the relationship being approximately 1.0, and an intercept near 0, the two devices can be judged as yielding similar results (Moore, et al., 1997). In the case of the current evaluation, the slope of the relationship between the ParvoMedics and PFS measurements of VO_2 did not significantly differ from 1.0, nor did the intercept differ from zero.

Interestingly, no differences were seen in expired pulmonary ventilation (V_E), even though the systems compute pulmonary ventilation using differing methods. Both the PFS and ParvoMedics use a differential pressure pneumotach to measure flow; however, the PFS measures flow on the inspiratory side, while the ParvoMedics measures flow on the expiratory side. The software of both systems uses the Haldane Transformation correction (Luft, 1973) to correct for the differences in inspiratory and expiratory ventilation when calculating VO_2 and VCO_2 . The PFS measures inspired ventilation, rather than the more traditional expired measurement, to minimize the chance of condensation or saliva contamination on the pneumotach surface. This design feature is important for spaceflight because traditional methods to reduce contamination from moisture contained in expired gas rely on gravity and thus would not work in the ISS environment.

Use of the PFS during exercise testing onboard the ISS should provide improved accuracy of the estimates of aerobic fitness which currently rely on the assumption that VO_2 per given workload is similar on ISS to preflight measurements. However, the prediction of peak oxygen uptake from the HR response to submaximal exercise, even with gas exchange measurements, will still have limited accuracy. If

NASA desires accurate aerobic capacity measurements, serious consideration must be given to the performance of maximal exercise testing with metabolic gas analysis.

In conclusion, laboratory evaluation of the PFS demonstrated that it provides similar results to those measured by the reference metabolic gas analysis system. It is recommended that the PFS be incorporated into the standard periodic fitness evaluation testing performed onboard the ISS.

5.0 REFERENCES

- Astrand, P.O. and Ryhming, I. A nomogram for calculation of aerobic capacity (physical fitness) from pulse rate during submaximal work. *J. Appl. Physiol.* 7:218, 1954.
- Bland, J.M. and Altman, D.G. Statistical methods for assessing agreement between two methods of clinical measurement. *Lancet I*: 307-310, 1986.
- Bassett, D.R., Howley, E.T., Thompson, D.L., King, G.A., Strath, S.J., McLaughlin, J.E., Parr, B.B. Validity of inspiratory and expiratory methods of measuring gas exchange with a computerized system. *J. Appl. Physiol.* 91:218-224, 2001.
- Carter, H., Jones, A.M., Maxwell, N.S., Doust, J.H. The effect of interdiurnal and diurnal variation on oxygen uptake kinetics during treadmill running. *J. Sports Sci.* 20: 901-9, 2002.
- Convertino, V.A. Exercise as a countermeasure for physiological adaptation to prolonged spaceflight. *Med. Sci. Sport. Exerc.* 28: 999-1014, 1996.
- Crouter, S.E., Antczak, A., Hudak, J.R., Della Valle, D.M., Haas, J.D. Accuracy and reliability of the ParvoMedics TrueOne 2400 and MedGraphics VO2000 metabolic systems. *Eur. J. Appl. Physiol.* 98: 139-151, 2006.
- Jensen, K., Jorgensen, S., Johansen, L. A metabolic cart for measurement of oxygen uptake during human exercise using inspiratory flow rate. *Eur. J. Appl. Physiol.* 87: 202-206, 2002.
- Lee, S. M. C., Moore Jr, A.D., Barrows, L.H. Variability of Prediction of Maximal Oxygen Consumption on the Cycle Ergometer Using Standard Equations. Houston, TX: Johnson Space Center, 1993, pp 1-8. (NASA TM-3412).
- Lepers, R., Millet, G.Y., Maffiuletti, N.A., Hausswirth, C., Brisswalter, J. Effect of pedaling rates on physiological response during endurance cycling. *Eur. J. Appl. Physiol.* 85: 392-395, 2001.
- Luft, U.C., Myrthe, L.G., Loeppky, J.A. Validity of Haldane calculation for estimating respiratory gas exchange. *J. Appl. Physiol.* 34: 864-865, 1973.
- Michel, E.L., Rummel, J.A., Sawin, C.F., Buderer, M.C., Lem, J.D. Results of Skylab medical experiment M171-metabolic activity. In: *Biomedical Results From Skylab*, R.S. Johnston and L.F. Dietlein, eds. NASA-SP-377, 372-387, 1977.
- Moore, A.D., Lee, S.M.C., Greenisen, M.C., Bishop, P. Validity of a heart rate monitor during work in the laboratory and on the Space Shuttle. *Am. Indust. Hygiene Assoc. J.* 58: 299-301, 1997.
- Shephard, R.J. Tests of maximum oxygen uptake: a critical review. *Sports Med.* 1: 99-124, 1984.
- Versteeg, P.G.A. and Kippersluis, G.J. Automated systems for measurement of oxygen uptake during exercise testing. *Int. J. Sports Med.* 10: 107-112, 1989.

REPORT DOCUMENTATION PAGE			Form Approved OMB No. 0704-0188	
Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project (0704-0188), Washington, DC 20503.				
1. AGENCY USE ONLY (Leave Blank)	2. REPORT DATE September 2007	3. REPORT TYPE AND DATES COVERED NASA Technical Paper		
4. TITLE AND SUBTITLE Validation of the Pulmonary Function System for Use on the International Space Station			5. FUNDING NUMBERS	
6. AUTHOR(S) Frank A. McCleary, M.S., Alan D. Moore, Jr., Ph.D., Wyle Life Sciences, Houston, TX R. Donald Hagan., Ph.D., NASA Johnson Space Center, Houston, TX				
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) Lyndon B. Johnson Space Center Houston, Texas 77058			8. PERFORMING ORGANIZATION REPORT NUMBERS S-1006	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES) National Aeronautics and Space Administration Washington, DC 20546-0001			10. SPONSORING/MONITORING AGENCY REPORT NUMBER TP-2007-214756	
11. SUPPLEMENTARY NOTES				
12a. DISTRIBUTION/AVAILABILITY STATEMENT Available from the NASA Center for AeroSpace Information (CASI) 7115 Standard Hanover, MD 21076-1320 Category: 52			12b. DISTRIBUTION CODE	
13. ABSTRACT (Maximum 200 words) Aerobic deconditioning occurs during long-duration spaceflight despite the use of exercise countermeasures. As a part of International Space Station (ISS) medical operations, periodic tests designed to estimate aerobic capacity are performed prior to, during, and after missions of greater than 30 days in duration. These tests track changes in aerobic fitness and determine the effectiveness of exercise countermeasures. The purpose of this investigation was to compare exercise metabolic gas analysis measurements (including oxygen consumption) obtained by the Pulmonary Function System (PFS) to those collected using a reference metabolic gas analysis system: the ParvoMedics TrueOne© 2400 system (ParvoMedics, Salt Lake City, UT). This system has been extensively validated and is currently utilized by the NASA's Exercise Physiology Laboratory for pre- and post-flight testing astronauts assigned to ISS flights. Laboratory evaluation of the PFS demonstrated that it provides similar results to those measured by the reference metabolic gas analysis system. It is recommended that the PFS be incorporated into the standard periodic fitness evaluation testing performed onboard the ISS.				
14. SUBJECT TERMS Exercise physiology, fitness, ergometers, pulmonary functions			15. NUMBER OF PAGES 22	16. PRICE CODE
17. SECURITY CLASSIFICATION OF REPORT Unclassified	18. SECURITY CLASSIFICATION OF THIS PAGE Unclassified	19. SECURITY CLASSIFICATION OF ABSTRACT Unclassified	20. LIMITATION OF ABSTRACT Unlimited	
